


# MOBILE ASPIRATORS

Models 3040, 4042

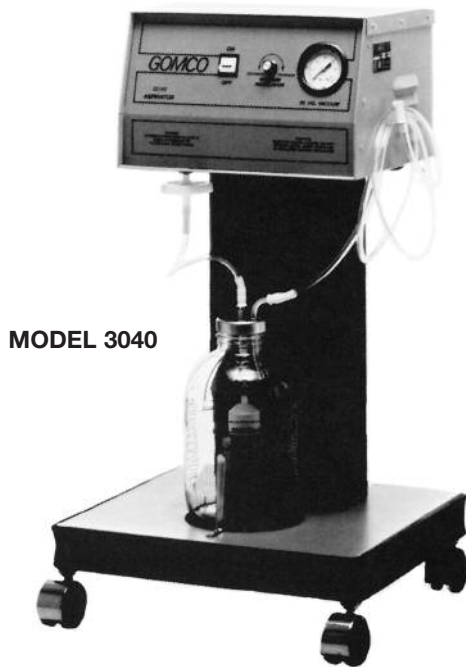

**MOBILE ASPIRATORS** – Gomco® stand-mounted units offer dependable, easy to regulate suction for hospital use. Stands include 3" non-marking casters, handy cord storage, chemical resistant stainless steel top, compact stand design, new float shut-off assembly, graduated one-gallon polycarbonate plastic collection bottle or 2100 ml disposable canister, bacteria filter to help prevent fluid and aerosol contamination, and a hospital grade plug.

**MODELS 3040/4042** – These units feature an insulated hidden motor and a durable stainless steel top. Model 3040 features a 1/8 horsepower diaphragm motor, and is UL listed and ULC certified. Model 4042 utilizes a 1/10 horsepower rotary compressor motor, and is UL listed.

Suction systems are equipped with an easy-to-read gauge, calibrated in both inches Hg and mm Hg, and a suction regulating valve. These models also include illuminated on/off switches.

MODEL	COLLECTION BOTTLE	MOTOR	COMPRESSOR	VACUUM RANGE	FLOW RATE	DIMENSIONS HxDxW	SHIPPING WEIGHT
3040	Polycarbonate Plastic	1/8 hp	Diaphragm	0-22 in. Hg	26 lpm	33 x 17 x 18 in.	70 lbs.
3041	Disposable			0-559 mm Hg		84 x 43 x 46 cm	32 kg
4042	Polycarbonate Plastic	1/10 hp	Carbon	0-25 in. Hg	30 lpm	33 x 17 x 18 in.	76 lbs.
4043	Disposable		Rotary Vane	0-635 mm Hg		84 x 43 x 46 cm	34 kg

## ORDERING INFORMATION

w/Polycarbonate Plastic Bottle		
	w/Disposable Canister	Voltage
01-21-3040	01-21-3041	110/50
01-22-3040	01-22-3041	115/60
01-62-3040*	01-62-3041*	115/60
01-63-3040	01-63-3041	220/50
01-64-3040	01-64-3041	220/60
01-12-4042	01-12-4043	110-115/50-60
01-52-4042*	01-52-4043*	115/60

\*Tested to 60601-1 IEC Standard